


Rapid Test Guide

HSV 1


Type: Herpes Simplex Virus 1
 Product: Rapid Test Kit
 Specimen: Whole Blood



Clean your finger. Use the alcohol pad to wipe the tip of your finger.

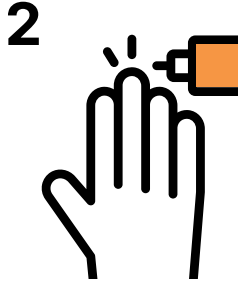


Prick your finger. Use the lancet to prick your fingertip and gently squeeze until a drop of blood appears.

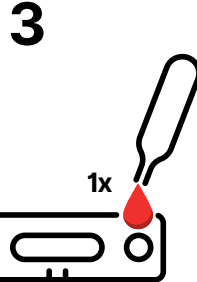


Add blood to the cassette. Put 1 drop of blood into the small round hole marked "S" on the test cassette.

Note: Make sure you add enough blood for the test to work. If the test strip doesn't start to get wet within 1 minute, add one more drop of buffer to the same spot.

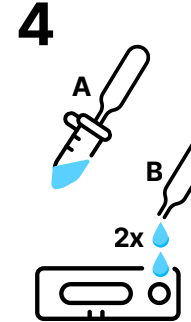


To add buffer liquid:

- Use pipette to get buffer from tube.
- Drop 2 drops into the round test well where blood was placed.
- Tube has extra, no need to use it all.

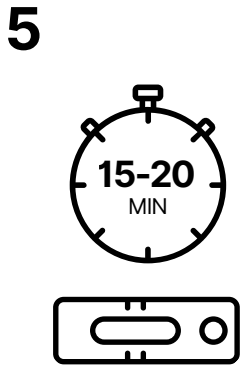

Wait for 15 minutes. Do not read the result after 20 minutes.

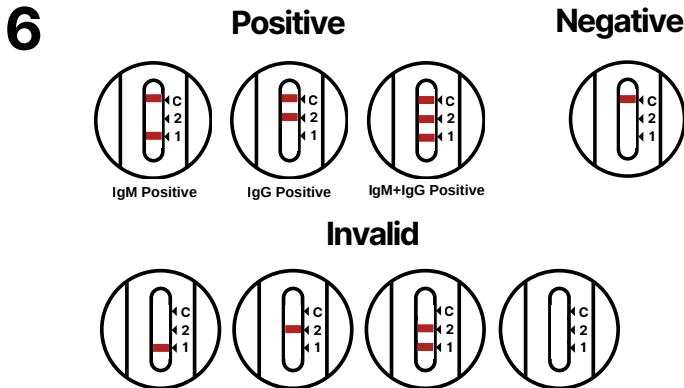

Check the result: 3 lines (C+2+1) or 2 lines (C+1 or C+2) = Positive. 1 line at C = Negative. No line at C = Invalid, test again.



Throw away the used kit properly. Put all used parts in a covered trash bin.


Rapid Test Guide

HSV 2


Type: Herpes Simplex Virus 2
 Product: Rapid Test Kit
 Specimen: Whole Blood



Clean your finger. Use the alcohol pad to wipe the tip of your finger.

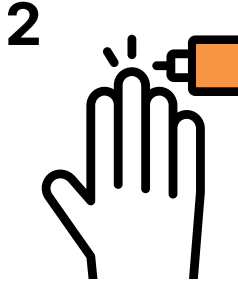


Prick your finger. Use the lancet to prick your fingertip and gently squeeze until a drop of blood appears.

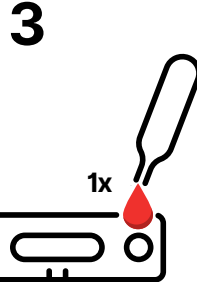

Add blood to the cassette. Put 1 drop of blood into the small round hole marked "S" on the test cassette.

Note: Make sure you add enough blood for the test to work. If the test strip doesn't start to get wet within 1 minute, add one more drop of buffer to the same spot.

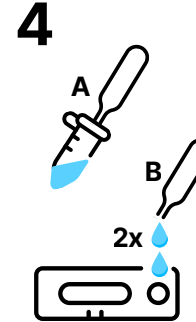

To add buffer liquid:

- Use pipette to get buffer from tube.
- Drop 2 drops into the round test well where blood was placed.
- Tube has extra, no need to use it all.

5



Wait for 15 minutes. Do not read the result after 20 minutes.

6

Positive



IgM Positive




IgG Positive






IgM+IgG Positive

Negative

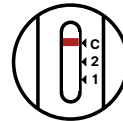

Invalid

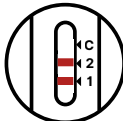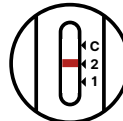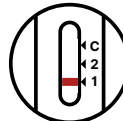

Check the result: 3 lines (C+2+1) or 2 lines (C+1 or C+2) = Positive. 1 line at C = Negative. No line at C = Invalid, test again.

7



Throw away the used kit properly. Put all used parts in a covered trash bin.